



Implant3D

Dental Implant Simulation

CT Scan Protocol

Patient ID

Name

Birth Date

Sex

M

F


1. Make CT Scan with tilt angle equal to zero.
2. Create CD-ROM with **Axial only** in **DICOM 3.0** format.
3. Separate Upper Jaw and Lower Jaw scans in two CD-ROM or two folders in same CD-ROM.

Best Size Pictures:

512x512 pixels

Best Axials thickness:

from 0.3 to 0.6 mm

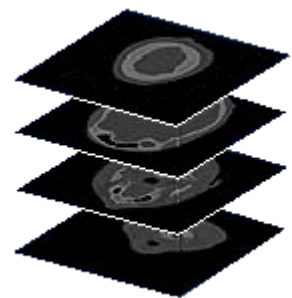

Please give this document to the Patient for Radiologist.